

Vita



Neonatal incubator





Intensive care of preterm or sick babies should not only secure the **survival**. Also devotion and comfort must be given to the baby in order to make **life** possible. The incubator **Vita** will help you to achieve this purpose. It conjoins the stability of a closed ambience with the good access and contact possibilities offered by open care systems.

Always a good climate

Air and skin temperature, humidity and oxygen concentration are selected according to the baby's need and are servo controlled. A brand-new air flow technique warms the patient evenly by minimizing the transepidermal humidity losses.

Ergonomic design


The perspex double-wall assembly has four ports. For quick access all walls can be folded down with one hand, the canopy can be removed easily. For passing tubes and cables sufficient grommets are provided.

For taking X-ray it is not necessary to open the wall assembly, the tray for the cassette is inserted from the outside. Of course, for special procedures the bed can be pulled-out completely.

The incubator is also available with an integrated weighing scale. There is an electromotive bed tilt in both directions with automatic stop in weighing position. The desired working height can be adjusted by foot pedal.

Storage space and choice of accessories

A large shelf and two drawers are provided for storing material necessary for treatment. The accessories from the THERMOCARE and VARIOTHERM devices are compatible and can be fixed at the rails and the spar assembly.

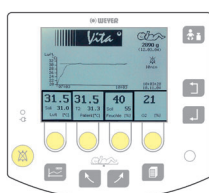

The control module

Innovative technology controls and monitors the selected values for air and patient temperature as well as for humidity and oxygen concentration. The selected parameters and measured values are indicated clearly structured with a graphic display of the development during the last three hours.

All values can be retraced and evaluated for the period of one week.

